

PRO 30

YOUR POWERFUL ALLROUNDER
FOR DENTAL LABS



rapidshape

PRO 30

Powerful Performance

- 4k Light Engine printing with $\pm 25\mu\text{m}$ pixel size optimized for high-precision printing
- Powerful printing speeds with patented Force Feedback and RS Turbo technology
- Built-in chamber heating allowing for resin flexibility

Scalable Printing for Dental Labs

200+ validated resins available for 15+ dental applications



Crown & Bridge Models
up to 6 Crown & Bridge Models

25 minutes

50 μm layer thickness



Denture Bases
up to 13 Denture Bases

90 minutes

100 μm layer thickness



Surgical Guides
up to 10 Surgical Guides

30 minutes

100 μm layer thickness

Reliable

- **Integrated Camera** reporting print job success to RS Customer Portal for production control and quality assessment
- Built to last DLP 3D printer with high-quality glass reservoirs
- Intelligent connectivity from software, to 3D Printer to WASH to CURE

Rapid Support

1-Year CARE Plus Plan included with purchase of PRO 30 3D printer

Service includes:

- 90-Day Coach with Rapid Shape Digital Workflow Specialist
- Continued Support and access to workflow guides and troubleshooting tools
- Annual sensor send-in
- Full repair coverage & immediate device replacement if necessary



A.RS.PRO30_RS

RS Crystal Polish

High-Shine Splints Requiring Minimal Polishing



RS Crystal Polish Reservoir
Item Number: BGM183401



Compatible Materials

detax keystone
industries

Up to 28 Splints printed in 90 min.

- Crystal clear print results with no yellowing effect after WASH & CURE cycle
- Reduces up to 90 % of polishing required after curing



PRO30 Accessories

Accessory Package

Item Number: RS006175

Your Setup for Continuous, Clean Printing

- Includes UV light proof storage for up to 5 reservoirs
- 4x Glass Reservoirs with covers
- Additional consumables



PRO 30 Reservoir

Item Number: BGM183400



Durable Glass Reservoir

- Tempered glass reservoirs without consumption limitation
- RFID chip integrated to indicate resin type

Post-Processing Solution

Fully streamlined connectivity between Rapid Shape 3D printers, WASH and CURE removing guesswork and manual parameter input from start to finish



WASH



Item Number: A.RS.PROWASH_RS2

- Max chamber capacity: 6 Full arch print jobs
- 2-step wash cycle and drying mechanism
- Reduced smell with activated carbon filtration
- Validated and pre-installed wash parameters for 200+ resins
- Stackable with CURE

CURE

Item Number: A.RS.PROCURE_RS2

- Max chamber capacity: 6 full arch print jobs
- Integrated vacuum system eliminating need for Nitrogen tanks or generators
- Achieving peak performance 3D prints with 360° curing in a controlled vacuum environment
- Nitrogen inlet built in

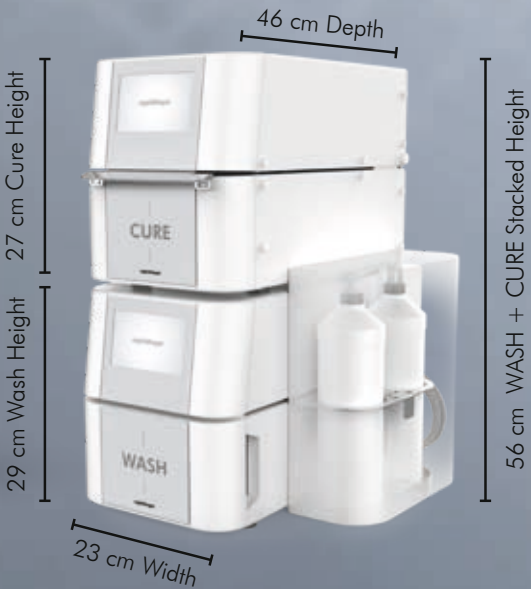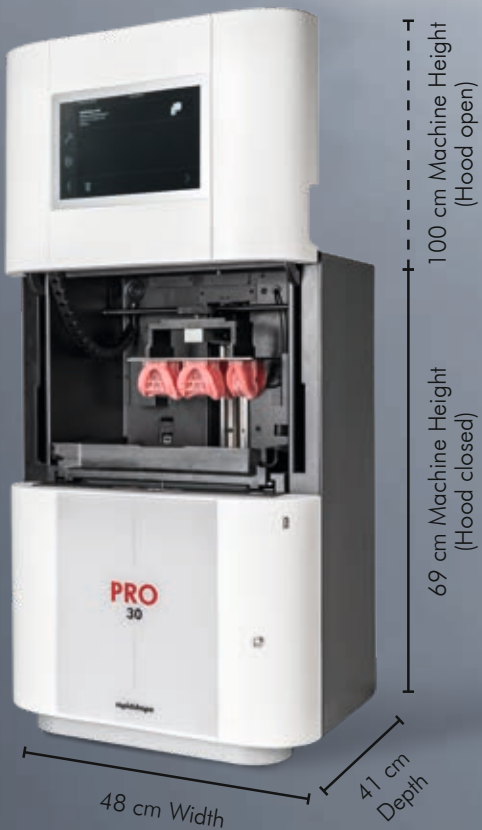

SCAN TO VIEW WASH



SCAN TO VIEW CURE



rs VIVO



RS VIVO Model

Premium model resin ideal for esthetic and precision dental models, including die and implant models.

Available colors: Ivory and Sand.

RS VIVO Model Natural Pink

High-quality 3D printing resin designed to visually replicate the natural appearance of oral soft tissue. Ideal for shade matching for veneers, crowns, and other prosthetic work.

Available color: Natural Pink



PRODUCT	SIZE	ITEM NUMBER
RS VIVO Model Ivory	1 kg	RS006190
RS VIVO Model Ivory	5 kg	RS006216
RS VIVO Model Sand	1 kg	RS006192
RS VIVO Model Natural Pink	1 kg	RS006191



printed with **rs** VIVO
Model Natural Pink and Ivory

Book A Demo



Scan to Schedule
remote 30-min demo
with RS Digital Workflow
Specialists

Contact Information for Customer Support:
service@rapidshape.de